

Quick Start QRS 101



1. Connect power cord and applicators to the top end of the control unit.



2. Turn on the power.



3. Select applicator A or B (large mat or small pad).



4. Select treatment program. i.e. Relax, Basis, Vital.



5. Select treatment time, i.e. 8, 16, or 24 minutes. Press the button for 2 seconds and use the arrow buttons to decrease or increase by 1 minute increments.



6. Select intensity, i.e. Sensitive, Medium, or Intensive. Press the button for 2 seconds and use the arrow buttons to decrease or increase the intensity.



7. Press Start. At the end of the treatment, the system will automatically turn off.



8. A light is attached to each applicator that will light up during treatment.

Quick Start QRS 101



1. Connect power cord and applicators to the top end of the control unit.



2. Turn on the power.




3. Select applicator A or B (large mat or small pad).



4. Select treatment program. i.e. Relax, Basis, Vital.



5. Select treatment time, i.e. 8, 16, or 24 minutes. Press the button for 2 seconds and use the arrow buttons to decrease or increase by 1 minute increments.

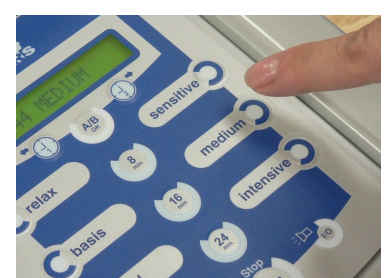

6. Select intensity, i.e. Sensitive, Medium, or Intensive. Press the button for 2 seconds and use the arrow buttons to decrease or increase the intensity.



7. Press Start. At the end of the treatment, the system will automatically turn off.



8. A light is attached to each applicator that will light up during treatment.

For more detailed information, please refer to the instruction manual.

For more detailed information, please refer to the instruction manual.